


Instruction for FeLV Ag/FIV Ab test



Sample fill line (black line) : 1 drop (10 μ l)

1

Sampling & loading

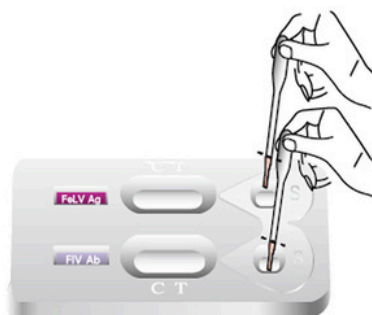

**X 2
(two times)**


2 drops of sample in case of FeLV Ag

1 drop of sample in case of FIV Ab

•• Do not overload the sample !!

2

Absorption

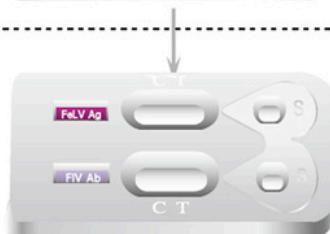


Wait for complete absorption

3

Buffer addition

00:00:00

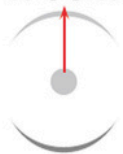

2 drops of buffer

•• after complete absorption




4

Reading

00:05:00~00:10:00



+

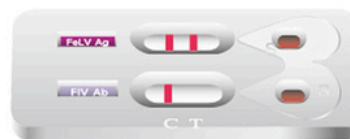

Read the result between 5 ~ 10 min

-



5

Discard

00:10:00



Discard the used device
In valid result after 10 min

