

Instruction for Giardia Ag test



The amount of sample : 100 μl


1

Sampling

1. Do not use sticky or a lump of feces.
2. Don't collect too much feces.

Maximum quantity : 30 mg

(Fig : Swab with 30 mg of feces)

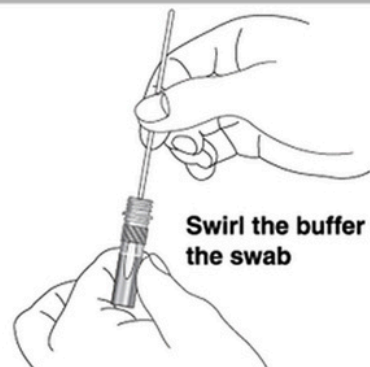

Swirl the buffer with the swab

2

Use the Supernatant layer



Wait until feces particles sink to the bottom



Supernatant layer for use

Feces particle layer



3

Loading

- 4 drops (100 μl) of the diluted sample
- **Don't overload the sample!!**



4

Reading

00:05:00~00:10:00



+



Read the result between 5~10 min

-



5

Discard

00:10:00



: Discard the used device

Invalid result after 10min

